

30" Ductless Fume Hood

Model # SS-330-DCH

MADE IN
THE USA



SINCE 1983

This **30" Wide Ductless Fume Hood** combines powerful air flow with a sturdy and compact enclosure for optimal respiratory and environmental protection. This fume hood protects the operator and environment from hazardous fumes or particulate generated from applications performed within the hood.

This **Ductless Fume Hood's** powerful fan directs fumes into the filter chamber, composed of a variety of filter media selected depending on the application. This filter media includes **HEPA** [up to 99.97% efficiency on particles 0.3 microns and larger], **ASHRAE** [up to 95% efficiency on particles 0.5 microns and larger], **Activated Carbon**, and **specialty-blended filtration** [*i.e. Acid Gas, Mercury, Aldehyde, Ammonia*].

Typical uses for this **Ductless Fume Hood** include *chemical fumes, pharmaceutical compounding, soldering applications, dust removal, biological applications, and a variety of applications requiring fume and dust removal*. The **Ductless Fume Hood** can provide an effective and economical solution for many commercial, medical, laboratory and industrial applications that require a small containment area. For additional visibility, a transparent version is available at extra cost.

Our line of ductless fume hoods are not intended for sterile compounding applications.



30" Ductless Fume Hood



30" DCH provides fume control for a various medical applications.

Product Features


- **Includes:** variable speed controller (115V only), velometer, and LED light
- Quiet operation
- Long filter life
- Simple, quick "no tool" filter change
- Portable and lightweight
- Reliable, low maintenance operation
- **Standard Sizes:** 12", 18", 24", 30", 40", 50", 60", 70"
- Custom Sizes Available
- **Optional:** clear anti-static vinyl curtains, Magnehelic Gauge, hour counter, and spill tray.



21221 FM 529 Rd, Cypress, TX 77433
1.800.799.4609 • sales@sentryair.com

www.sentryair.com

Product Specifications

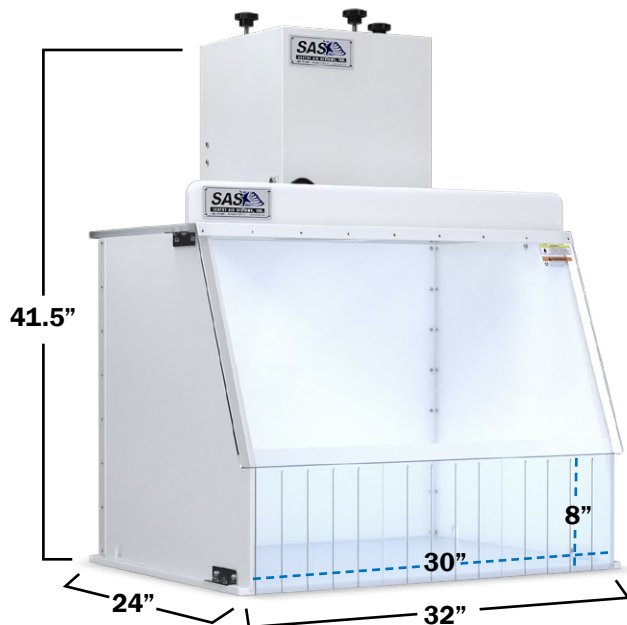

Shown with outer dimensions.
Dimensions are approximate.

Hood Dimensions (D x W x H)	Outer: 24" x 32" x 25" Interior: 23.5" x 30" x 24" Access Area: 8" x 30"
Hood Material	Sides: 0.25" white acrylic Base: 0.5" white high density polyethylene Top: 0.5" clear polycarbonate Front Window: 0.25" clear polycarbonate
Blower Dimensions	15.5" L x 13" W x 14.5" H
Blower Material	16 Ga. carbon steel with polyester epoxy coating on all surfaces
Weight (approximately)	100 lbs
Air Volume (approximately)	350 CFM <i>(varies with filter media)</i>
Inlet Velocities on High (approximately)	Carbon Pre-filter + HEPA Filter – 130 FPM high, 25 FPM low Carbon Pre-filter + Heavy Duty Carbon Filter – 115 FPM high, 20 FPM low
Electrical	115/1/160, 2 amps 220/1/50, 0.75 amps 8' grounded power cord with NEMA 5-15P plug
Sound Level	61 to 66 dba (low to high)
Maximum Temperature	Operating temperature not to exceed 170° F
Main Filter Options	HEPA (up to 99.97% efficient on particles down to 0.3 microns) ASHRAE (up to 95% efficient on particles down to 0.5 microns) Activated Carbon Specialty-blended filter media (i.e. acid gas, mercury, aldehyde, ammonia)
Warranty (From Date of Shipment)	Limited two-year warranty on defects due to materials or workmanship.
Patent #	5,843,197



21221 FM 529 Rd
Cypress, TX 77433
1.800.799.4609
sales@sentryair.com

www.sentryair.com